

# Training on Proficiency Testing Scheme GeneXpert DTS

## Module 4: Workflow on preparation of DTS panels

Venue  
Presenter  
Date

# Introduction

- To prepare a quality panels internal developed SOPs should be followed and adhere to in the preparation process.

# Objectives

By the end of this module participants should be able to:

- Describe the technical work flow process in preparation of genexpert PT panels
- Describe the management work process in the preparation of genexpert PT panels

# Module outline

- Workflow on preparation of DTS panels( Technical)
- Workflow on preparation of DTS panels(management)

# Work flow Processes (Technical)

- Culture Selection and Freezer Storage
- Inactivation and Stock Preparation
- Pre-testing and Stock Selection
- Panel Aliquoting
- Panel Validation
- Labeling and Packing
- Shipping



# 1. Culture selection and freezer storage

Isolates selected should have all the parameters tested by Xpert

- susceptible
- resistant
- Negative

And from the subculture isolates a freezer stock should be prepared for future use



# Prepare culture isolates

- Inoculate the selected freezer stocks on MGIT, Incubate until MGIT positive
- Remove positive MGIT, write the date of positivity on the side of tube and place in axillary incubator at 35-37°C for 4-6 days.
- Store at 2-8°C until inactivation

# Prepare culture isolates



Inoculate seed culture  
into MGIT tube



Incubate in MGIT  
until positive



Incubate 4-6 days  
longer



Store at 2-8° until  
inactivation





## 2. Inactivation


This process is undertaken to ensure non viable organisms reas used for the preparation of the panel

Different equipment can be used for this process ie


- Hot air oven 85° C( Preferred since the temperatures are closely controlled and monitored)
- Water bath 85° C
- Autoclave 121 Pa

The identified isolates on MGIT are used and after inactivation they can be pooled together

# Briefly vortex MGIT cultures to be inactivated

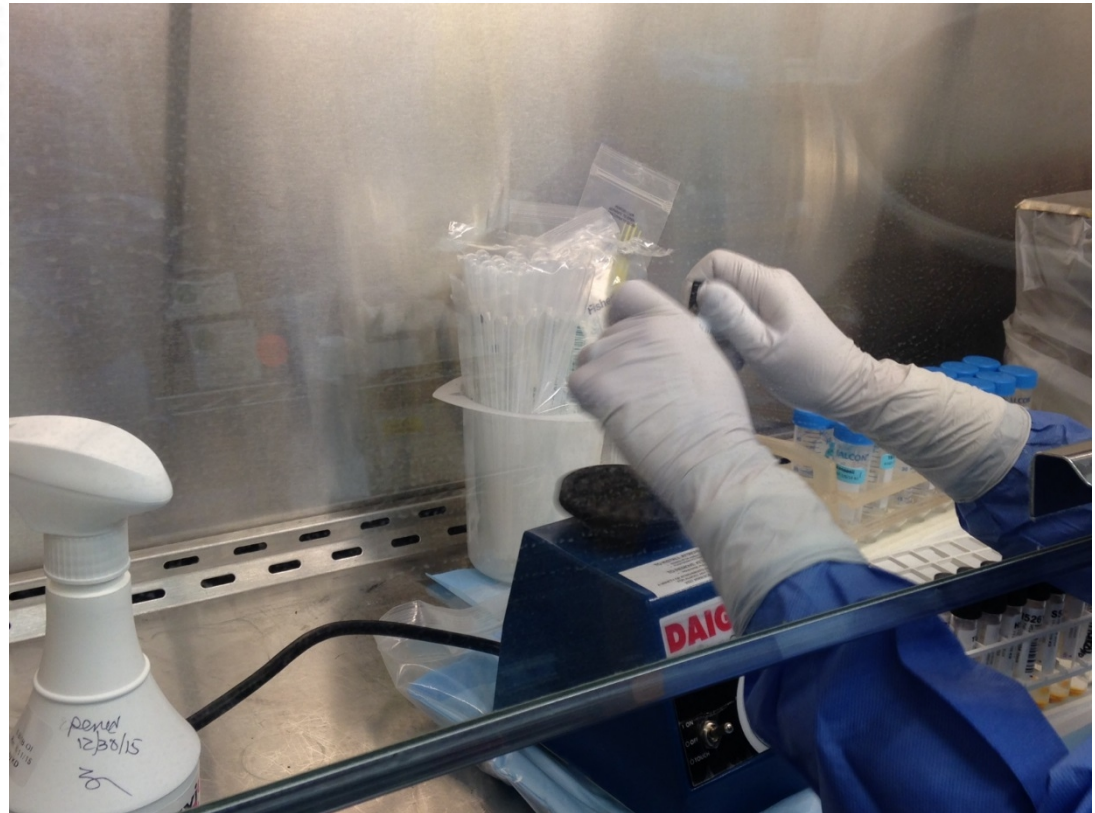

# Inactivation verification

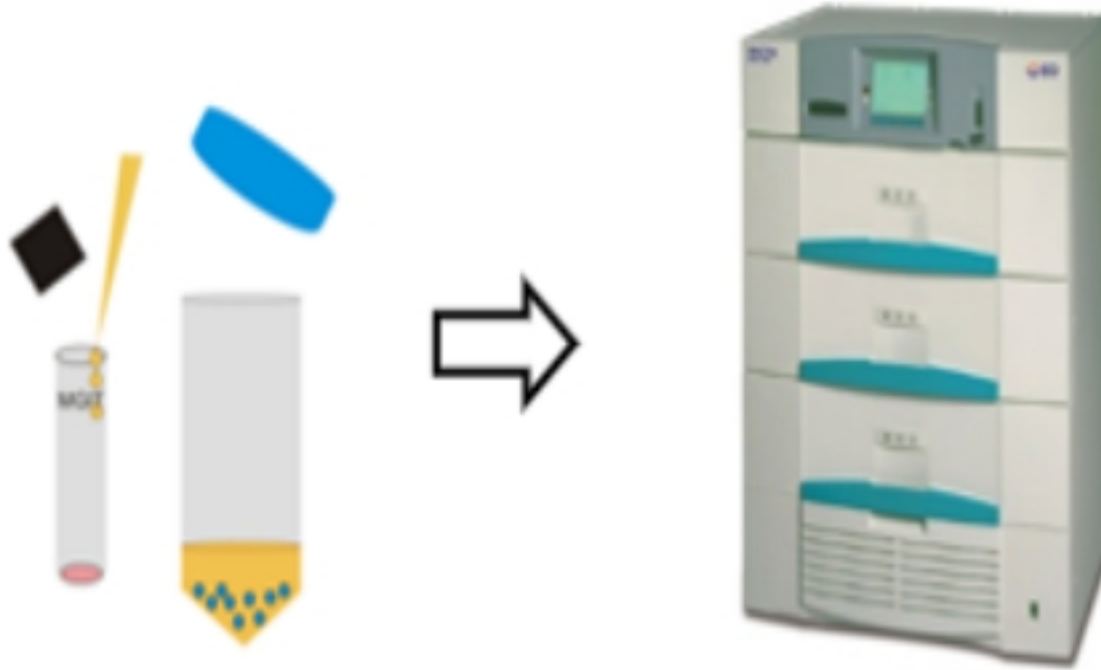

**Incubate stock solution in MGIT for 84 days**


# 3. Pretesting and stock selection

The inactivated organisms are pretested to see if

🌐 The dilution factor gives the expected results

🌐 Choose the best stocks based on the factor above.

# Dilute and aliquot

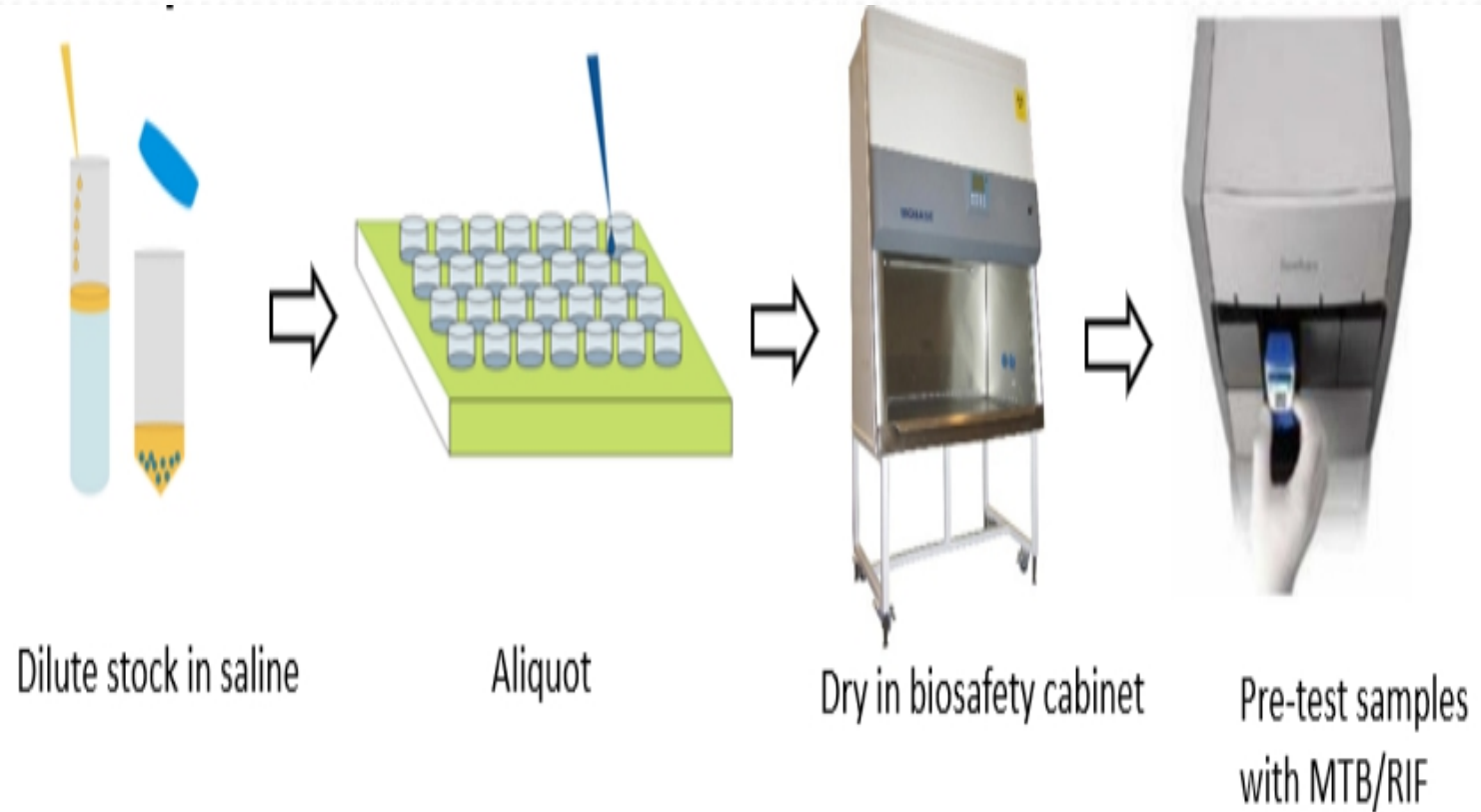



# 4. Panel aliquoting

Aliquoting the final panel depends on

- 🌐 The pretest results

- 🌐 Inactivation verification results( any growth on this process leads to automatic disqualification of the stock)



# 5. Panel validation

5% of the aliquoted stock is used for the validation process

Essential to ensure that the aliquoted stocks give a

 To obtain consensus results

 Taken as a test of homogeneity

 Stability of the panel

# 6. Labelling and packaging

Good quality labels should be considered

- To avoid falling off or fading and facilitate testing the actual panels and obtaining correct results

Quality packaging materials should be considered

- To maintain the integrity of the panels during transportation


# 7. Shipping culture PT panels

Local and international regulations should be followed to ensure PT panels are delivered to the testing facilities

🌐 In time

🌐 safely

# Dried Tube Specimens using Heat Inactivation

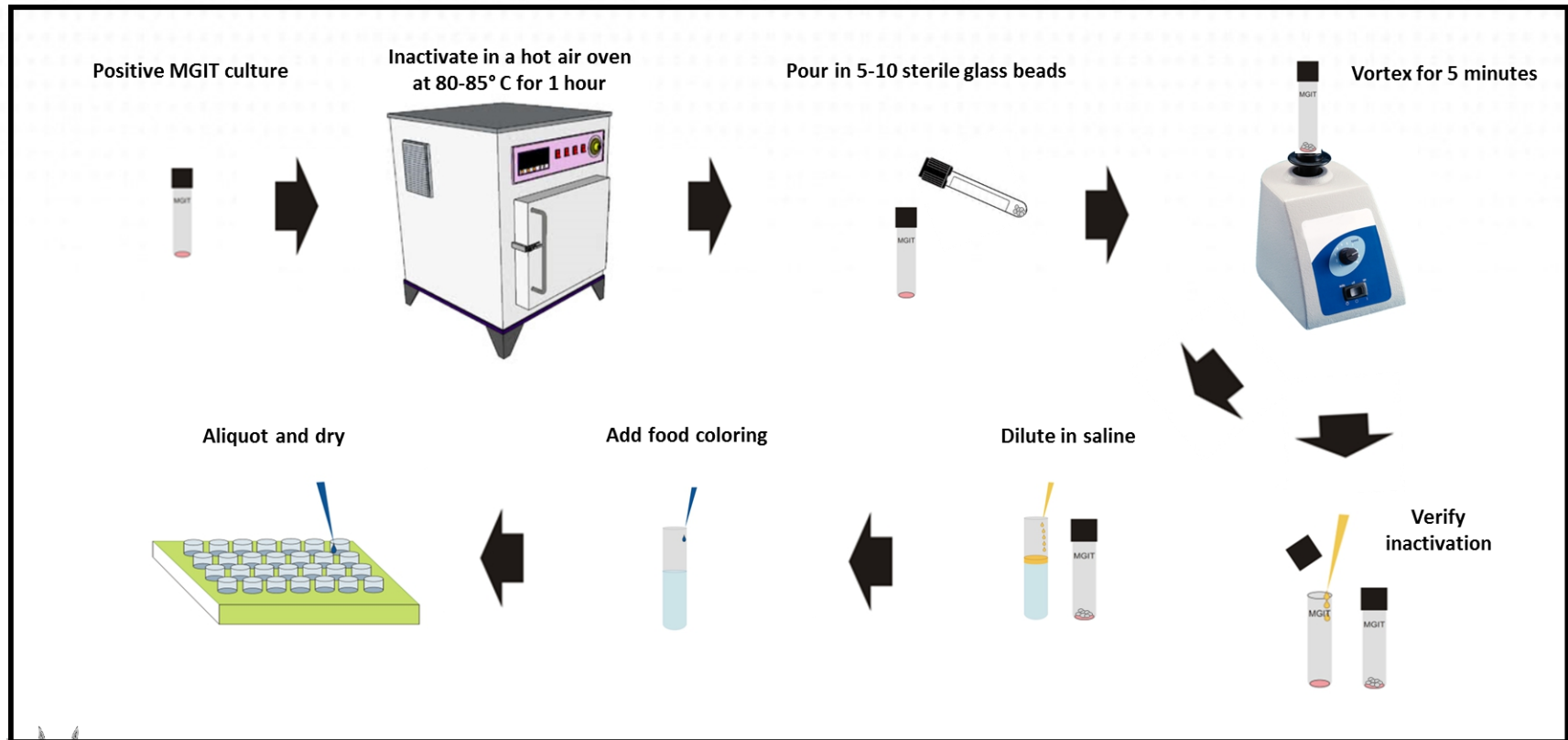

# Workflow processes( program management)

- **Program Scheduling and Enrollment**
- **Data Collection and Analysis**
  - Online and Manual (MS Excel and Word)
- **Following Up Non-Conformities with Sites**
  - Contacting sites and documenting root cause analysis
  - Corrective actions
- **Record keeping for genexpert PT scheme**
- **Indicators genexpert PT scheme**
- **Check ISO 17043:10 indicators where applicable**



# Assessment

1. List the 6 work flow technical processes in preparation of genexpert PT panels
2. What is the method employed to verify inactivation of the inoculum?
3. List atleast 4 programme/ scheme management workflow process in the genexpert PT scheme?



# Summary

- The technical work flow process include, Culture Selection Freezer Storage, Inactivation, and Stock Preparation, Pre-testing, Stock Selection, Panel Aliquoting, Panel Validation, Labeling and Packing
- MGIT culture is the method employed to verify inactivation of the inoculum
- MGIT culture is preferred to LJ for growth of the inoculum because of the growth-time factor



# REFERENCES

- ISO 13528:2005, *Statistical methods for use in proficiency testing by interlaboratory comparisons*
- ISO Guide 34, *General requirements for the competence of reference material producers*
- ISO Guide 35, *Reference materials – General and statistical principles for certification*
- Guide 34, ISO Guide 35 and ISO 13528 (homogeneity and stability)
- ISO/IEC Guide 98-3, *Uncertainty of measurement – Part 3: Guide to the expression of uncertainty in measurement* (GUM:1995)
- ISO/IEC 17011:2004, *Conformity assessment – General requirements for accreditation bodies accrediting conformity assessment bodies*
- ISO/IEC 17025, *General requirements for the competence of testing and calibration laboratories*



# Acknowledgments

