FACILITATOR GUIDE

THE TB SAMPLE REFERAL SYSTEM (TSRS) TRAINING

MODULE 6

Module 6 TB sample packaging and transportation: Concept of triple packaging Sample transportation Waste disposal.


Acknowledgements

This Module was prepared by Akum Mary.

Version 1.0 Effective date: 01-Jun-2019

TSRS/FG/006 Page **1** of **4**

Part 6.0: **TB sample packaging and transportation:** Concept of triple packaging Sample transportation Waste disposal.

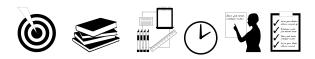

OBJECTIVE: To understand the concept of triple packaging, labelling of sample boxes during shipment and waste disposal.

MATERIALS:

- 1. Handouts,
- 2. Slides
- 3. Computer
- 4. Over head projector
- 5. Flip chart
- 6. Markers and pens
- 7. Note books
- 8. Materials for Practical session on completing the laboratory request forms

TIMELINE: 45 Minutes

METHODOLOGY:

- 1. Lecture
- 2. Discussion
- 3. Practical session on completing the laboratory request forms

ADVANCE PREPARATION:

- 1. Printing of handouts,
- 2. Familiarize oneself with sample labeling and completing lab request form

FACILITATORS STEP-BY STEP INSTRUCTIONS:

- 1. Welcome and Introduction
- 2. Present module overview
- 3. Ask question about module expectations
- 4. Start the practical presentation
- 5. Allow participants to carry out Practical session and observe for to ensure correct
- 6. Ask if there is any question.
- 7. Distribute handouts

FACILITATORS NOTES

- 1. Take critical understanding of triple packaging, UN labels recommended for labeling of biological material and waste disposal.
- 2. Engage the trainees frequently by giving NTRL examples and ask them how they do it in their facilities too.

SLIDE OF POWERPOINT PRESENTATION

TSRS/FG/006 Version 1.0

Effective date: 01-Jun-2019



SITUATION ANALYSIS/ EXERCISES

Break class into groups of 3-5 people

1. How would you package a sputum sample to be transported from laboratory A to B for Xpert MTB/RIF assay?

ASSESSMENT REVIEW

- 1. What do you understand by triple packaging?
- 2. How do you handle wastes generated during sample packaging
- 3. What is the standard TAT in days for transporting samples between the peripheral lab and the central lab?

REFERENCES

- GLI TB training package http://www.stoptb.org/wg/gli/trainingpackages.asp
- Global tuberculosis report 2019, https://apps.who.int/iris/bitstream/handle/10665/329368/9789241565714-eng.pdf?ua=1
- The Handbook Laboratory Diagnosis of Tuberculosis by Sputum Microscopy http://www.stoptb.org/wg/gli/assets/documents/TB%20MICROSCOPY%20HA
 http://www.stoptb.org/wg/gli/assets/documents/TB%20MICROSCOPY%20HA
 https://www.stoptb.org/wg/gli/assets/documents/TB%20MICROSCOPY%20HA
 https://www.stoptb.org/wg/gli/assets/documents/
 https://www.stoptb.org/wg/gli/assets/
 https://www.stoptb.org/wg/gli/assets/
 https://www.stoptb.org/wg/gli/assets/
 https://www.stoptb.org/wg/gli/assets/
 https://www.stoptb.org/wg/gli/
- TUBERCULOSIS LABORATORY BIOSAFETY MANUAL https://apps.who.int/iris/bitstream/handle/10665/77949/9789241504638_e
 ng.pdf?sequence=1

Version 1.0 Effective date: 01-Jun-2019

TSRS/FG/006 Page 4 of 4